



Patient Care Furniture



ISO CERTIFIED 13485:2016 & 9001:2015



Indtech expertise in manufacturing innovative, comfortable and easy-to-use healthcare furnitures with utmost safety for patients in mind. Working for frontline healthcare workers including doctors and nurses and patients who seeks care and comfort in their difficult times, Indtech never compromise quality.

Our products ranges from ICU advanced care and ICU support to Dialysis, Examination support, Ward care, General support, Transfer support and Mother care. The right talent in production applied upfront allows us to meet almost all of our client's requirements in time.

QUALITY POLICY

The quality policy of IndTech Health Care Pvt. Ltd., is to provide products that meet or exceed customer expectations, comply with all applicable regulations, continually improve in everything we do, and maintain the effectiveness of the quality management system.

WHAT MAKES US SPECIAL



CUSTOMISED PRODUCTS

We provide our customers with the products that meet their specific demands and requirements.



ISO 13485:2016 CERTIFIED

The products are accredited to ISO 13485:2016 & 9001: 2015 and CE certification for European Economic Area



INTERNATIONAL STANDARDS

Every time we make sure that the products we manufacture meets the international export standards.



PRODUCT MAINTENANCE

We don't quit once the product is handed over to our clients. We provide maintenance services for the products on site.



ON TIME DELIVERY

We make sure that the products we manufacture get delivered to our clients on time.



STATE-OF-THE-ART PRODUCTION SYSTEM

We use the latest technologies to manufacture the products and ensure the utmost quality.

In the face of laborious working environment in hospitals, Indtech offers comfortable and easy to handle furnitures with utmost safety for patients. With varieties of different designs and grades from deluxe to normal, we meet requirements for each and every hospitals. Check out our products to know more.





LEARN MORE ABOUT INDTECH PATIENT CARE FURNITURES.

ICU Advanced Care
Dialysis
Ward Care
Transfer Support
Mother Care

ICU Support
Examination Support
General Support



IHC 1000 - FIVE FUNCTION FULL MOTORIZED ICU BED

- **Headboard & Footboard (ABS):** Easily detachable for easy intubations and resuscitation.
- **Electric operation:** Provided for all five functions that can be operated with a hand held remote.
- **Side railings:** Split and swing-down type made of PP plastic materials.
- **Emergency CPR position** assisted by manual operation.
- **Nursing control panel** at the bows
- **Bed top:** Four sections made of CRCA steels with perforations for sufficient airflow.
- **SS Saline rod & Provisions** on all corners.
- **Buffers:** Provided on four corners to the bed from damages on collisions and movements.
- **Urine bag holder** provision at suitable location.
- **Base & Frame:** Made of CRCA tubes.
- **Wheels:** 125mm non marking castor with central locking system.

BED DIMENSIONS

Length : 2100mm
Width : 1060mm

HEIGHT ADJUSTMENT


Min : 505mm without mattress.
Max : 750mm without mattress.

MATTRESS PLATFORM DIMENSIONS

Length : 2030 mm
Width : 910 mm

ADJUSTABLE ANGLE RANGE (DEGREE)


Backrest angle : 0 - 85
Knee rest angle : 0 - 40
Trendelenburg angle : 0 - 17
Reverse Trendelenburg angle : 0 - 17

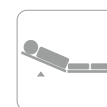


Backrest




Knee rest

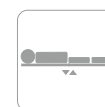

Hi-Low

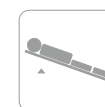

Trendelenburg

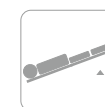

R. Trendelenburg



CPR zero positioning



IHC 1001 - FIVE FUNCTION FULL MOTORIZED ICU BED

- **Headboard & Footboard (ABS):** Easily detachable for easy intubations and resuscitation.
- **Side railings:** Split and swing-down type made of PP plastic materials.
- **Electric operation:** Provided for all five functions that can be operated with a hand held remote.
- **Bed top:** Four sections made of CRCA steel sheets with perforations for sufficient airflow.
- **SS Saline rod & Provisions on all corners.**
- **PU Buffers:** Provided on four corners to the bed from damages on collisions and movements.
- **Urine bag holder** provision at suitable location.
- **Pre-treated and Epoxy powder coat finish.**
- **Base and frame** made of CRCA tubes.
- **Wheels:** Four 125mm non marking castors with diagonally opposite individual locking system.

BED DIMENSIONS

Length : 2100mm
Width : 1060mm

HEIGHT ADJUSTMENT

Min : 505mm without mattress.
Max : 750mm without mattress.

MATTRESS PLATFORM DIMENSIONS


Length : 2030 mm
Width : 910 mm

ADJUSTABLE ANGLE RANGE (DEGREE)

Backrest angle : 0 - 85
Knee rest angle : 0 - 40
Trendelenburg angle : 0 - 17
Reverse Trendelenburg angle : 0 - 17

OPTIONAL ACCESSORIES

Nurse control system
Mattress
Reading / under bead light
Central lock system

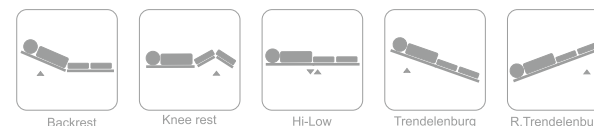



IHC 1002 - FIVE FUNCTION FULL MOTORIZED ICU BED

- **Headboard & Footboard (ABS):** Easily detachable for easy intubations and resuscitation.
- **Side railings:** Split and swing-down type made of PP plastic materials.
- **Emergency CPR position** assisted by manual operation.
- **Nursing control panel** at the bows
- **Electric operation:** Provided for all five functions that can be operated with a hand held remote.
- **Bed top:** Four sections made of CRCA steel tube, Removable mattress platform of suitable size made up of HDPE top
- **SS Saline rod & Provisions** on all corners.
- **PU Buffers:** Provided on four corners to the bed from damages on collisions and movements.
- **Urine bag holder** provision at suitable location.
- **Pre-treated and Epoxy powder coat** finish.
- **Base and frame** made of CRCA tubes.
- **Wheels:** Four 125mm non marking castors with central locking system.

BED DIMENSIONS

Length : 2100mm
Width : 1060mm

HEIGHT ADJUSTMENT

Min : 485mm without mattress.
Max : 686mm without mattress.

MATTRESS PLATFORM DIMENSIONS

Length : 2000 mm
Width : 900 mm

ADJUSTABLE ANGLE RANGE (DEGREE)

Backrest angle : 0 - 85
Knee rest angle : 0 - 40
Trendelenburg angle : 0 - 17
Reverse Trendelenburg angle : 0 - 17

OPTIONAL ACCESSORIES


Nurse control system
Mattress
Reading / under bead light
Central lock system




Backrest




Knee rest

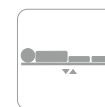

Hi-Low

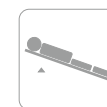

Trendelenburg

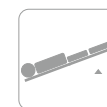

R. Trendelenburg



CPR zero positioning



IHC1005 - FIVE FUNCTION MANUAL ICU BED

- **Headboard & Footboard (ABS):** Easily detachable for easy insulation and resuscitation.
- **Side railings:** Split and swing-down type made of PP plastic materials.
- **Manual operation:** Regulated with high quality cranks for head-rest, knee-rest, hi-low, trendelenburg & reverse trendelenburg positions.
- **Bed top:** Four sections made of CRCA steel sheets with perforations for sufficient airflow.
- **SS Saline rod & Provisions** on all corners.
- **PU Buffers:** Provided on four corners to the bed from damages on collisions and movements.
- **Urine bag holder** provision at suitable location.
- **Pre-treated and Epoxy powder coat** finish.
- **Base and frame** made of CRCA tubes.
- **Wheels:** Four 125mm non marking castors with diagonally opposite individual locking system.

BED DIMENSIONS

Length : 2100mm
Width : 1060mm

HEIGHT ADJUSTMENT

Min : 550mm without mattress.
Max : 810mm without mattress.

MATTRESS PLATFORM DIMENSIONS

Length : 2030 mm
Width : 910 mm

ADJUSTABLE ANGLE RANGE (DEGREE)

Backrest angle : 0 - 85
Knee rest angle : 0 - 40
Trendelenburg angle : 0 - 17
Reverse Trendelenburg angle : 0 - 17

OPTIONAL ACCESSORIES

Mattress
Central lock system





IHC 1006 - THREE FUNCTION FULL MOTORIZED ICU BED

- **Headboard & Footboard (ABS):** Easily detachable for easy intubations and resuscitation.
- **Side railings:** Collapsible type made of Aluminium.
- **Electric operation:** Provided for all three functions that can be operated with a hand held remote.
- **Bed top:** Four sections made of CRCA steel sheets with perforations for sufficient airflow.
- **SS Saline rod & Provisions** on all corners.
- **Urine bag holder** provision at suitable location.
- **Pre-treated and Epoxy powder coat** finish.
- **Base and frame** made of CRCA tubes.
- **Wheels:** Four 125mm non marking castors with diagonally opposite individual locking system.

BED DIMENSIONS

Length : 2100mm
Width : 1060mm

ADJUSTABLE ANGLE RANGE (DEGREE)


Backrest angle : 0 - 85
Knee rest angle : 0 - 40

HEIGHT ADJUSTMENT

Min : 505mm without mattress.
Max : 750mm without mattress.

MATTRESS PLATFORM DIMENSIONS

Length : 2030 mm
Width : 910 mm

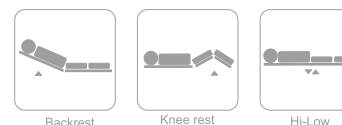



IHC1007 - THREE FUNCTION MANUAL ICU BED

- **Headboard & Footboard (ABS):** Easily detachable for easy insulation and resuscitation.
- **Side railings:** Collapsible type made of Aluminium.
- **Manual operation:** Regulated with high quality cranks for head-rest, knee-rest, hi-low positions.
- **Bed top:** Four sections made of CRCA steel sheets with perforations for sufficient airflow.
- **SS Saline rod & Provisions** on all corners.
- **Urine bag holder** provision at suitable location.
- **Pre-treated and Epoxy powder coat** finish.
- **Base and frame** made of CRCA tubes.
- **Wheels:** Four 125mm non marking castors with diagonally opposite individual locking system.

BED DIMENSIONS

Length : 2100mm
Width : 1060mm

ADJUSTABLE ANGLE RANGE (DEGREE)


Backrest angle : 0 - 85
Knee rest angle : 0 - 40

HEIGHT ADJUSTMENT

Min : 550mm without mattress.
Max : 810mm without mattress.

MATTRESS PLATFORM DIMENSIONS

Length : 2030 mm
Width : 910 mm

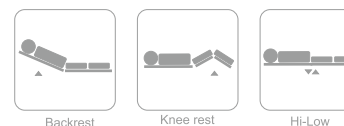



IHC1008 - THREE FUNCTION MANUAL ICU BED (MOVEMENT COT)

- **Headboard & Footboard (ABS):** Easily detachable for easy intubations and resuscitation.
- **Side railings:** Collapsible type made of Aluminium.
- **Manual operation:** Regulated with high quality cranks for back-rest, knee-rest, trendelenburg and reverse trendelenburg elevations.
- **Bed top:** Four sections made of CRCA steel sheets with perforations for sufficient airflow.
- **SS Saline rod & Provisions** on all corners.
- **Urine bag holder** provision at suitable location.
- **Pre-treated and Epoxy powder coat** finish.
- **Base and frame** made of CRCA tubes.
- **Wheels:** Four 125mm non marking castors with diagonally opposite individual locking system.

BED DIMENSIONS

Length	: 2100mm
Width	: 1060mm
Height	: 610mm

MATTRESS PLATFORM DIMENSIONS


Length	: 1850 mm
Width	: 850 mm

ADJUSTABLE ANGLE RANGE (DEGREE)

Backrest angle	: 0 - 85
Knee rest angle	: 0 - 40
Trendelenburg angle	: 0 - 17
Reverse Trendelenburg angle	: 0 - 17

OPTIONAL ACCESSORIES

Mattress

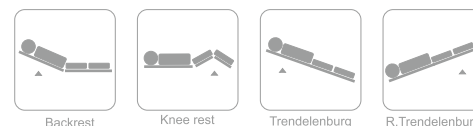



IHC 1016 - RECOVERY BED

- **Headboard & Footboard (ABS):** Easily detachable for easy intubations and resuscitation.
- **Manual operation:** Regulated with high quality cranks for back-rest, trendelenburg and reverse trendelenburg elevations.
- **Bed top:** Two sections made of CRCA steel sheets with perforations for sufficient airflow.
- **SS Saline rod & Provisions** on all corners.
- **Urine bag holder** provision at suitable location.
- **Pre-treated and Epoxy powder coat** finish.
- **Base and frame** made of CRCA tubes.
- **Legs:** Mounted on high quality rubber bushes.

BED DIMENSIONS

Length : 2100mm
Width : 1060mm
Height : 610mm

MATTRESS PLATFORM DIMENSIONS


Length : 1850 mm
Width : 850 mm

ADJUSTABLE ANGLE RANGE (DEGREE)


Backrest angle : 0 - 85
Knee rest angle : 0 - 40

OPTIONAL ACCESSORIES


Mattress
Urine bag holder

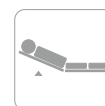

Backrest

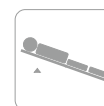

Trendelenburg

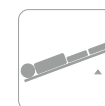

R.Trendelenburg



IHC 1050 - DIALYSIS CHAIR

- CRCA pipe frame work with three section top.
- Electric operation: Back-rest, Foot-rest, hi-low and reverse trendelenburg positions with four actuators.
- Adjustable arm-rest.
- Adjustable foot-rest.
- SS Saline rod.
- Pre-treated and Epoxy powder coated.
- Wheels: Four 125mm non marking castors with diagonally opposite individual locking system.



DIMENSIONS

Length : 1925mm
Width : 925mm

ADJUSTABLE ANGLE RANGE (DEGREE)

Backrest angle : 0 - 80
Knee rest angle : 0 - 40

HEIGHT ADJUSTMENT


Min : 650mm.
Max : 780mm.

OPTIONAL ACCESSORIES

Central lock system

UPHOLSTERY DIMENSIONS

Length : 2030 mm
Width : 910 mm

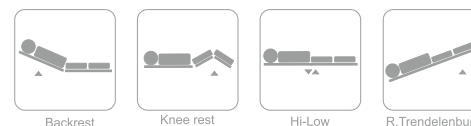



IHC1101 - MOTORIZED FOWLER BED

- **Headboard & Footboard (ABS):** Easily detachable for easy intubations and resuscitation.
- **Side railings:** Collapsible type made of Aluminium.
- **Electric operation:** Provided for back-rest and knee-rest functions that can be operated with a hand held remote.
- **Bed top:** Four sections made of CRCA steel sheets with perforations for sufficient airflow.
- **SS Saline rod & Provisions** on all corners.
- **Pre-treated and Epoxy powder coat finish.**
- **Base and frame** made of CRCA tubes.
- **Wheels:** Four 125mm non marking castors with diagonally opposite individual locking system.

BED DIMENSIONS

Length : 2100mm
Width : 1060mm
Height : 610mm

MATTRESS PLATFORM DIMENSIONS


Length : 1850 mm
Width : 850 mm

ADJUSTABLE ANGLE RANGE (DEGREE)


Backrest angle : 0 - 85
Knee rest angle : 0 - 40

OPTIONAL ACCESSORIES

Mattress
Urine bag holder

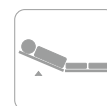

Backrest

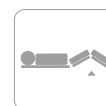

Knee rest



IHC1102 - FOWLER BED

- **Headboard & Footboard (ABS):** Easily detachable for easy intubations and resuscitation.
- **Side railings:** Split and swing-down type made of PP plastic materials.
- **Manual operation:** Regulated with high quality cranks for back-rest and knee-rest elevations.
- **Bed top:** Four sections made of CRCA steel sheets with perforations for sufficient airflow.
- **SS Saline rod & Provisions** on all corners.
- **Pre-treated and Epoxy powder coat** finish.
- **Base and frame** made of CRCA tubes.
- **Wheels:** Four 125mm non marking castors with diagonally opposite individual locking system.

BED DIMENSIONS

Length : 2100mm
Width : 1060mm
Height : 610mm

MATTRESS PLATFORM DIMENSIONS


Length : 1850 mm
Width : 850 mm

ADJUSTABLE ANGLE RANGE (DEGREE)


Backrest angle : 0 - 85
Knee rest angle : 0 - 40

OPTIONAL ACCESSORIES

Mattress
Urine bag holder

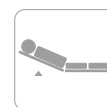

Backrest

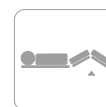

Knee rest



IHC1103 - FOWLER BED

- **Headboard & Footboard (ABS):** Easily detachable for easy intubations and resuscitation.
- **Side railings:** Collapsible type made of Aluminium.
- **Manual operation:** Regulated with high quality cranks for back-rest and knee-rest elevations.
- **Bed top:** Four sections made of CRCA steel sheets with perforations for sufficient airflow.
- **Saline rod provisions** on four corners.
- **Pre-treated and Epoxy powder coat finish.**
- **Base and H-frame** made of CRCA tubes.
- **Wheels:** Four 125mm non marking castors with diagonally opposite individual locking system.

BED DIMENSIONS

Length : 2100mm
Width : 1060mm
Height : 610mm

MATTRESS PLATFORM DIMENSIONS


Length : 1850 mm
Width : 850 mm

ADJUSTABLE ANGLE RANGE (DEGREE)


Backrest angle : 0 - 85
Knee rest angle : 0 - 40

OPTIONAL ACCESSORIES

Mattress
Urine bag holder
Telescopic SS Saline rod

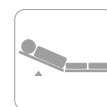

Backrest

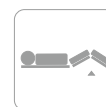

Knee rest




IHC1105 - FOWLER BED SS

- **Head & Foot** made of stainless steel tubes.
- **Side railings:** Swing-away type made of stainless steels.
- **Manual operation:** Regulated with high quality cranks for back-rest and knee-rest elevations.
- **Base and H-frame** made of stainless steel tubes.
- Available in 304 and 202 grades.
- **Legs:** Mounted on high quality rubber bushes.
- **Saline rod provisions** on four corners.

MATTRESS PLATFORM DIMENSIONS

Length : 1854mm
Width : 864mm
Height : 584mm

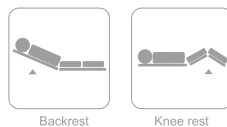


IHC1106 - FOWLER BED MS

- **Head & Foot** made of steel tubes.
- **Side railings:** Swing-away type made of MS.
- **Manual operation:** Regulated with high quality foldable cranks for back-rest and knee-rest elevations.
- **Pre-treated and Epoxy powder coated.**
- **Base and H-frame** made of CRCA steel tubes.
- **Legs:** Mounted on high quality rubber bushes.
- **Saline rod provisions** on four corners.

MATTRESS PLATFORM DIMENSIONS

Length : 1854mm
Width : 864mm
Height : 584mm

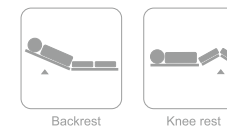



IHC 1110 –ELECTRIC SEMI FOWLER BED

- **Headboard & Footboard (ABS):** Easily detachable.
- **Bed top** made with ABS material
- **Side railings:** Split and swing-down type made of PP plastic materials.
- **Electric operation:** Provided for back-rest function that can be operated with a hand held remote.
- **Pre-treated and Epoxy powder coat** finish.
- **Base and frame** made of CRCA steel tubes.
- **Wheels:** Four 125mm non marking castors with diagonally opposite individual locking system.

MATTRESS PLATFORM DIMENSIONS

Length : 1854mm
Width : 864mm
Height : 584mm



IHC 1111 – SEMI FOWLER BED

- **Headboard & Footboard (ABS):** Easily detachable.
- **Side railings:** Collapsible type made of Aluminium.
- **Manual operation:** Regulated with high quality cranks for back-rest elevation.
- **Pre-treated and Epoxy powder coat** finish.
- **Base and H-frame** made of CRCA steel tubes.
- **Wheels:** Four 125mm non marking castors with diagonally opposite individual locking system.

MATTRESS PLATFORM DIMENSIONS

Length : 1854mm
Width : 864mm
Height : 584mm






IHC 1112 – SEMI FOWLER BED

- **Head & Foot** made of stainless steel tubes.
- **Side railings:** Swing-away type made of stainless steels.
- **Manual operation:** Regulated with high quality for back-rest elevations.
- **Base and H-frame** made of stainless steel tubes.
- Available in 304 and 202 grades.
- **Legs:** Mounted on high quality rubber bushes.

MATTRESS PLATFORM DIMENSIONS

Length : 1854mm
Width : 864mm
Height : 584mm

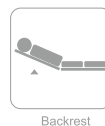

IHC 1113 – SEMI FOWLER BED

- **Head & Foot** made of MS.
- **Side railings:** Swing-away type made of MS.
- **Manual operation:** Regulated with high quality cranks for back-rest elevations.
- **Pre-treated and Epoxy powder coat finish.**
- **Base and H-frame** made of CRCA steel tubes.
- **Legs:** Mounted on high quality rubber bushes.

MATTRESS PLATFORM DIMENSIONS

Length : 1854mm
Width : 864mm
Height : 584mm





IHC 1114 – SEMI FOWLER BED

- **Head & Foot** made of stainless frame & wooden boards.
- **Side railings:** Swing-away type made of stainless steels.
- **Manual operation:** Regulated with high quality cranks for back-rest elevations.
- **Pre-treated and Epoxy powder coat** finish.
- **Base and H-frame** made of CRCA steel tubes.
- **Legs:** Mounted on high quality rubber bushes.

MATTRESS PLATFORM DIMENSIONS

Length : 1854mm
Width : 864mm
Height : 584mm



IHC 1115 – SEMI FOWLER BED

- **Head & Foot** made of MS.
- **Side railings:** Swing-away type made of MS.
- **Manual operation:** Regulated with high quality cranks for back-rest elevations.
- Removable commode
- **Pre-treated and Epoxy powder coat** finish.
- **Base and H-frame** made of CRCA steel tubes.
- **Legs:** Mounted on high quality rubber bushes.

MATTRESS PLATFORM DIMENSIONS

Length : 1854mm
Width : 864mm
Height : 584mm





IHC 1119 – PEDIATRIC BED

- **Bed top:** Perforated CRCA sheets.
- **Side rails:** Drop down type in full length.
- **Frame:** Rectangular CRCA tubes of approx 50mm x 25mm.
- **Pre-treated and Epoxy powder coat finish.**
- **Legs:** Mounted on high quality rubber bushes.



IHC 1120 – NORMAL HOSPITAL BED

- **H-Frame and Top** made of CRCA tubes.
- **Pre-treated and Epoxy powder coat finish.**
- **Legs:** Mounted on high quality rubber bushes.

Length : 1828mm

Width : 914mm

Height : 584mm

MATTRESS PLATFORM DIMENSIONS

Length : 1550mm

Width : 765mm

Height : 584mm



IHC 1822 – ATTENDANT BED

- **H-Frame and Top** made of CRCA tubes.
- **Pre-treated and Epoxy powder coat finish.**
- **Legs:** Mounted on high quality rubber bushes.

Length : 1830mm

Width : 610mm



IHC 1823 - ATTENDANT BED CUM CHAIR

- Manual adjustable back and leg-rest controlled by ratchet.
- Mechanical lever for back and leg rest operations
- Pre-treated and epoxy powder coated.
- Manual Chair position.
- Manual Relax position.



TECHNICAL SPECIFICATIONS

Length : 1850mm

Width : 520mm

Height : 545 mm

Backrest angle : 0°, 40°, 70°

Knee rest angle : 0°, 45°, 90°

Upholstery thickness : 75 mm

Safe working load : 130 Kg